



University of HUDDERSFIELD

University of Huddersfield Repository

Jones, Benjamin and Anson, T

In diamond health: diamond-like coatings for medical devices

Original Citation

Jones, Benjamin and Anson, T (2008) In diamond health: diamond-like coatings for medical devices. *Materials World*, 16 (8). pp. 24-26. ISSN 0967-8638

This version is available at <http://eprints.hud.ac.uk/id/eprint/19312/>

The University Repository is a digital collection of the research output of the University, available on Open Access. Copyright and Moral Rights for the items on this site are retained by the individual author and/or other copyright owners. Users may access full items free of charge; copies of full text items generally can be reproduced, displayed or performed and given to third parties in any format or medium for personal research or study, educational or not-for-profit purposes without prior permission or charge, provided:

- The authors, title and full bibliographic details is credited in any copy;
- A hyperlink and/or URL is included for the original metadata page; and
- The content is not changed in any way.

For more information, including our policy and submission procedure, please contact the Repository Team at: E.mailbox@hud.ac.uk.

<http://eprints.hud.ac.uk/>

In Diamond Health

Diamond-like coatings for medical devices

BJ Jones^{1}, T. Anson^{1,2}*

*1. Experimental Techniques Centre, Brunel University, Uxbridge, Middlesex UB8 3PH
UK*

2. Diameter Ltd., Brunel University, Uxbridge, Middlesex UB8 3PH UK

In contrast to the sparkling gemstone, diamond-like carbon (DLC) is a thin film coating that is dense, inert, low friction and hard wearing. The diverse range of areas exploiting DLC, including electronics, razors, manufacturing and medical devices, use a spectrum of materials that the name covers.

The umbrella term includes amorphous carbon, hydrogenated and tetrahedral variants – materials that can be likened to polymers, graphite or diamond. However, in contrast to the well defined crystalline structures of more well known carbon-based materials, DLC has no long-range order, but contains both graphitic and diamond-like bonded carbon atoms, and usually has a significant level of hydrogen. Tuning the balance of these constituents can finely tailor the properties of a material to match the application.

Interdisciplinary research, involving materials scientists, physicists, mechanical engineers, biomedical specialists and clinicians, is continuously expanding the potential and applications of DLC and enhanced carbon-based materials in the medical sector.

Diffusion barriers

Medical implants may contain undesirable elements in bulk or trace quantities. Nickel (Ni), chromium (Cr), cobalt (Co) and vanadium (V) can be toxic or cause sensitivity in some patients and are contained in commonly used metal alloys such as stainless steel (Fe-Ni-Cr), Ti-Al-V and Ni-Ti shape-memory alloy. Generally, at least one of the constituent metals is reactive in the presence of oxygen and readily forms protective oxide layers on the implant. These enable large numbers of devices to be implanted without adverse effects on biocompatibility. These oxide layers are, in most cases, effective in reducing contact between tissue and medical implant. However, a significant percentage of the population are sensitive to nickel, for example, and ions are still able to diffuse through the oxide barrier, with potentially adverse reactions.

A diffusion barrier that has little or no effect on human tissue and reduces undesirable metal ion migration is therefore desirable. Amorphous DLC, plasma deposited as a uniform thin film coating, shows encouraging results as a diffusion barrier. Coating

of devices external to the body, such as DLC-coated metal watch straps, provides a barrier between the body and the allergen, eliminating reaction to the nickel content of the metal.

Coatings can even enhance the aesthetics, providing a 'sparkling black' finish (see images below), and can be applied to spectacle frames and other metal items that come in to contact with the skin. Unlike varnish-type or polymer products, they are less susceptible to wear that would otherwise re-expose the Ni-containing substrate and trigger a dermatological complaint.



Titanium-alloy total knee prosthesis, coated (above) and uncoated (Below). Images courtesy of Diameter Ltd, Stanmore Implants Worldwide, RNOH

Potentially more of a hazard is that a patient may be sensitive to a metal in a medical implant, such as an orthopaedic or cardiovascular device. Orthopaedic prostheses such as hips, knees and shoulders are comprised of a substantial amount of metal with a large surface area exposed to surrounding hard and soft tissue. Any sensitivity is often detected prior to an operation by giving patients a metal allergy test to circumvent an untoward response to the implant. Diamond-like carbon has been successfully applied to a number of implanted devices, including total knee prostheses for people sensitive to constituents of a titanium alloy.


In addition to protecting the patient from medical devices, the reverse is often required to guard against the harsh environment of the human body. For example, implanted electronic sensors that interface with the neurological system, such as cochlear implants, can react to changes in the environment, fluid contamination or humidity, but need to retain a level of contact with the body. A thin film coating that prevents ingress of fluids can help the device remain functional.

More generally, coatings can reduce the effects of body fluids on implants, which could otherwise lead to increased corrosion and accelerated structural failure. Diamond-like carbon protected implants have the potential to reduce revision operations to replace damaged devices.

On the surface

The popularity of DLC is not solely related to diffusion barrier properties.

The surface topological, tribological and chemical properties make the material beneficial in many medical areas. Cardiovascular applications exploit the haemocompatibility and anti-thrombogenic properties of DLC – the surface structure reduces the likelihood of blood clots when a material is inserted into the blood stream. Capitalising on this, DLC-coated arterial stents are commercially available, preventing closure of arteries and ensuring unobstructed blood flow.

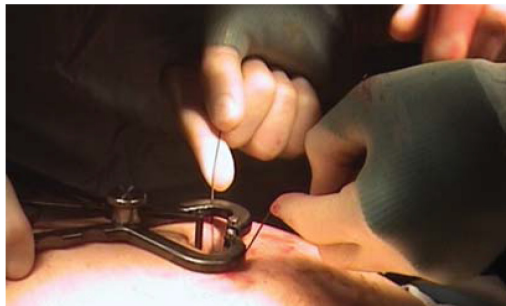

Top: Diamond-like carbon (DLC) coated aluminium based instrument used during gall-bladder removal. Image courtesy of Arun Mahendran, Hillingdon Hospital. Below: Partially DLC coated osteotomy screws. Image courtesy of Joe Franks, Diameter Ltd, RNOH

Urinary catheterisation is a common procedure for both short- and long-term use. Coating catheters with DLC reduces friction, both on insertion and in use, reducing local tissue damage and improving patient comfort. The catheter is an effective bacterial

pathway into the body and the coating helps to mitigate this, prohibiting biofilm formation and bacterial adhesion.

This property has led to further use of DLC in medical devices, such as in partially-coated osteotomy screws (see image left, bottom) used after the 7 July London bombings in 2005. These are transcutaneous prostheses that protrude outside the body and are used for small fracture fixation, such as hand repair after trauma. The DLC provides an inhospitable environment for bacterial colonisation and improves the rate of healing.

As well as remedial application to existing devices, DLC is facilitating new products which have the coating designed-in. A DLC-coated assist tool for gall bladder operations, developed at Brunel University, allows the use of aluminium for surgical tools, with the advantage of lighter weight and lower machining costs than traditional stainless steel. The DLC coating of the instrument dispels fears concerning metal transfer into the body.

Enhancements

In addition to fine tuning DLC and optimising it for the desired application, additional elements or structure can be added during or after deposition. Incorporating silicon into the films expands the range of materials that can be coated, improves adhesion and enhances haemocompatibility, as well as allowing increased flexing without fracture damage. This technique is used in a range of applications from drill bits to expanding stents.

Further processing techniques, such as fluorine incorporation during deposition, or post-deposition oxygen plasma treatment, can control surface hydrophobicity of the coating – a property that affects biocompatibility and fluid flow in medical devices, as well as water resistance of watches.



DLC coated watch straps. Images courtesy of Diameter Ltd

Work currently under development includes incorporation of biocides, such as macro-patterned or nano-particulate silver, into DLC. This dual coating combines the advantageous surface and diffusion barrier properties of DLC with the well-established biocidal effects of silver, and has the potential for use in surgical tools and biomedical implants to reduce the prevalence of hospital acquired infections.

Acknowledgements

Thanks to J. Franks and A. Mahendran.

Further Information

Dr Benjamin Jones,

Email: b.j.jones@physics.org